



# AC4 BED

Single | King Single

## USER MANUAL



Thank you for choosing this Aidacare AC4 Bed.  
This bed has been specifically designed to offer  
comfort, support and safety.

**It is important that you read this user manual  
carefully before using your new bed.**



This product features the CE mark, in compliance with the Medical Devices Regulation 2002 (as amended) Class 1, and the Medical Device Regulation 2017/745.



This product Features the UKCA mark, in compliance with Part II UK MDR 2002 (as amended) Class 1



**IMPORTANT:** Do not use the Aidacare AC4 Bed without having first read this user manual

© Aidacare Limited 2024

No part of this user manual may be copied, distributed, reproduced, scanned or stored on any electronic database, whether in whole or in part in any form or by any means.

Due care has been taken to ensure all the information contained in this user manual is correct at the time of printing. All measurements, pictures, colours and weight capacities are to be used as a guide only. We reserve the right to modify the design or appearance of any product displayed in this user manual without prior notice.

## TABLE OF CONTENTS

---

1. Symbols Used In This Manual	4
2. Symbols Used On The Product	5
3. Introduction	6
4. Safety Precautions	7
5. Electronic Functions	10
6. General Bed Operation	12
7. Specifications	19
8. Cleaning	20
9. Maintenance	21
10.Transport & Storage	23
11. Accessories	24
12.Warranty & Service	25

## 1. SYMBOLS USED IN THIS MANUAL

The symbols below are used throughout this owner's manual and on the product to identify warnings and important information. It is very important for you to read them and understand them completely.




Indicates a potentially hazardous condition/situation. Failure to follow designated procedures can cause either personal injury, component damage or malfunction. On the product, this icon is represented as a black symbol on a white triangle with a black border.




**ALWAYS!** These actions should be performed as specified. Failure to perform mandatory actions can cause personal injury and/or equipment damage. On the product, this icon is represented as a white infinity symbol on a black dot with a white border.



**DO NOT!** These actions are prohibited. These actions should not be performed at any time or in any circumstances. Performing a prohibited action can cause personal injury and/or equipment damage. On the product, this icon is represented as a white symbol with a black circle and black slash.

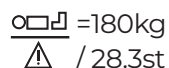

**MAXIMUM PATIENT WEIGHT.** This weight represents the total maximum weight of a patient on the bed. This weight (180kg / 28.3st) should not be exceeded at any time.

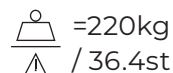


**SAFE WORKING LOAD (SWL).** This weight represents the total maximum weight that the bed can carry, including patient, mattress, bedding, bed ends and accessories. This weight (220kg / 36.4st) should not be exceeded at any time.



Please scan the QR code or visit our website **aidacare.com.au** and search 'Aidacare AC4 (BEB050400 or BEB050410)' to access product marketing and technical documentation.

## 2. SYMBOLS USED ON THE PRODUCT

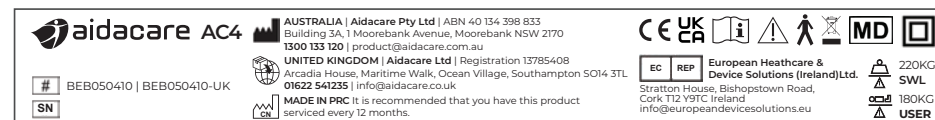

Sticker shown is not to scale.



Product Code



Serial Number



Manufacturer



Importer



Authorised Representative in EU



Country of Manufacture



UKCA Mark



Includes Electronics components and should be mixed with General household waste



CE Mark



Read User Manual



Caution



Medical Device



Safe Working Load



Maximum Patient Weight



Type B as per IEC60601-1



Class 2 Electrical Equipment

## 3. INTRODUCTION

---

The Aidacare AC4 Bed is an electrically operated Care Bed. It has a range of functions and positional movements that assist caregivers and provide for patient comfort. Ease of use, functionality, serviceability and cleaning practicality are key features of the bed. The bed has a fabricated steel structure with a powder coated finish. The plastic components are made from urethane, PP and POM nylon and ABS plastic.

### CLINICAL APPLICATIONS AND AREA OF USE

The bed is appropriate for adult patient care under supervised medical monitoring.

### KEY FEATURES

- Maximum safe working load (SWL) of the bed is 220kg (36.4st)
- Electrically adjusted height range of 300mm to 750mm
- 125mm diameter twin castors with central locking brakes and steering
- Four section electrically operated mattress platform with independent adjustment of backrest angle (maximum 70°) and knee section angle (maximum 40°)
- Mattress retainer handles
- Integrated foot end extension. Extends the bed by 195mm
- Clamp on  $\frac{3}{4}$  length folding side rails (optional accessory)
- Removable bed ends
- Wall protection system

## 4. SAFETY PRECAUTIONS

---

It is **important** that this User Manual is read in its entirety prior to operating the bed. The manual is designed to ensure that users are equipped with the information required to competently and responsibly operate the bed and its accessories.

It is recommended that the User Manual be kept in an easily accessible place for quick reference. If required Aidacare product specialists can provide additional training.

All warnings and cautions within this manual should be strictly adhered to. Failure to do so could result in harm to the patient or operator.

**Warnings** highlighted throughout this manual identify possible hazards in procedures or conditions which if not followed correctly could result in serious injury, or fatality.

**Cautions** highlighted throughout this manual identify possible hazards in procedures or conditions, which if not followed correctly could result in equipment damage, or failure.

### GENERAL

- Electrical equipment can be hazardous if used incorrectly.
- Residual current devices (RCD) are not supplied with the AC4 beds. Please consult with your Biomedical Engineer/advisor concerning RCD requirements.
- The clinically qualified person responsible should consider the age, size and condition of the patient before using the bed and its accessories.
- The maximum safe working load of the AC4 beds of 220kg (36.4st) includes the mass of the patient, mattress and any accessories fitted to the bed.
- Maximum mattress thickness that can be used on the AC4 beds is 170mm.
- Side rails must only be used with a mattress of appropriate size, which is approved for use with the bed.
- If a patient's condition could lead to a possible entrapment the mattress support platform should be left in the flat position whilst unattended.



**ALWAYS** position the bed at its lowest height when unattended by caregivers to minimise the risk of patient injury from falls when getting in and out of the bed



**ALWAYS** use the ExitSafe Button whenever the patient is entering or exiting the bed to ensure transfers are completed at a safe height to minimize falls risk



**ALWAYS** use a mattress of appropriate size and type. A range of suitable pressure reducing and pressure relieving mattresses are also available from Aidacare



**ALWAYS** make sure that the patient is safely positioned to avoid entrapment before operating the bed



**ALWAYS** ensure that attendant and patient limbs are free from potential entrapment areas before operating any functions



**ALWAYS** use the ExitSafe Button when treating a patient or transporting the bed to avoid unnecessary bending/stooping by care staff



**ALWAYS** apply the brakes when the bed is stationary



**ALWAYS** store the handset when not in use to avoid accidental operation of the beds functions. The handset can be stored by hooking it on to the side rail of the bed or positioning it in the handset holder bracket



**ALWAYS** ensure that all drips and drains (hospital lines) and power cables are free from possible entanglement with mechanical parts of the bed



**ALWAYS** transport/relocate beds using the EasyMove Transport System when uninhabited. This enables lift free transportation and reduces manual handling risks to staff



**DO NOT** place the handset, or power cord across the mattress surface of the bed when the bed is not connected to mains power. Always store the power cord using the fastening straps provided



**DO NOT** use electrically operated beds in the presence of flammable gases



**DO NOT** use the bed to support more than one patient. AC4 Beds are intended for single person use and may become unstable if used incorrectly. Usage outside of intended design will void all warranties associated with this product



**DO NOT** use accessories that are not designed, or approved for use with AC4 Beds

## 5. ELECTRONIC FUNCTIONS




**WARNING**

- Ensure the patient, caregivers and/or equipment is positioned to avoid being trapped or crushed during bed adjustment
- Ensure the power lead is not stretched excessively
- Ensure the handset and electrical leads cannot become entangled with moving parts of the bed
- DISCONNECT the power lead from the mains and stow before transporting the bed
- If the power lead shows any sign of damage, replace immediately
- Ensure safety key is removed when handset is not in use

### HANDSET BUTTON OPERATION

AC4 beds are powered by three electric actuators, which are controlled by the buttons on the handset. Adjustment occurs only while a button is pressed. Adjustment will stop when the button is released, or the movement selection reaches the end of its range of adjustment.

**CAUTION:** Do not operate the bed controls for more than the prescribed duty cycle of the control system. (See rate load page 18)

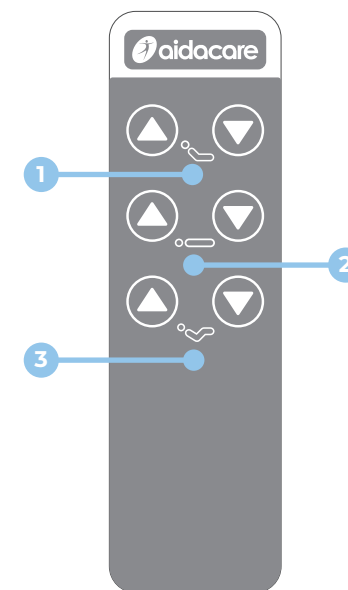

#### BACKREST ADJUSTMENT

1

To move the backrest up or down, press and hold the appropriate button on the handset. The backrest can be adjusted through a range of between 0° and 70°.

**WARNING:** Care should be taken whilst operating the backrest function to ensure hands and fingers are free from potential squeeze hazards.

#### MATTRESS PLATFORM HEIGHT ADJUSTMENT

2

To move the mattress platform up or down, press and hold the appropriate button on the handset. The platform height can be adjusted from a low of 300mm to a high of 750mm.

#### KNEE POSITIONING ADJUSTMENT

3

To move the knee positioning section up and down, press and hold the appropriate button on the handset. The knee section is adjustable through a range of 0° to 40°.

## 6. GENERAL BED OPERATION

### 6.1 CASTOR SYSTEM – BRAKES & STEERING

#### CENTRAL LOCKING CASTORS

AC4 Series beds have a central locking castor system that is operated with a full width foot pedal.

- With the pedal in the down **BRAKE** position, the brakes are applied to three castors of the bed simultaneously and the bed is immobilised.
- With the pedals in the horizontal **FREE** position, the brakes are released and the bed can be moved in any direction.
- With the pedals in the raised **STEER/TRACK** position, one steering castor at the head end of the bed is engaged to provide better control when pushing the bed in a straight line.

#### CENTRAL LOCK PEDAL POSITIONS



#### CAUTION:

- Disconnect the bed power cord before moving the bed.
- Do not remove the power cord from the source by pulling on the cord.
- Do not use the beds side rails to manoeuvre the bed.
- Do not pull on the power cord, or handset cord to manoeuvre the bed.
- Always apply the brakes when the bed is in a stationary position.

### 6.2 BED EXTENSION

The integrated bed extension can add an extra 200mm of length to the mattress platform.

#### MOVING BED EXTENSION FROM CLOSED/STANDARD BED POSITION

1. The bed extension is held in position by two threaded hand wheels on either side of the bed. Unscrew (anti clockwise) the hand wheels on both sides. It is not necessary to remove the hand wheels completely, unscrewing them by approximately 10mm of thread length will allow for the internal lock nut to travel freely past the thread of the hand wheel.
2. With both hand wheels now unscrewed, gently pull with an even force on both sides of the bed end. (2) Indicates the best position for easier extension if used on both sides of the panel at the same time.
3. Lock hand wheels in place on either side to lock bed end in.
4. Return the hand wheels back to their locked position (1)

**To reposition the bed end back to the closed/standard bed position use the same technique as above.**



**WARNING:** The bed extension is designed for use as a bed extension only and is not meant to support the body weight of a person sitting. The extra length of the bed with extra load can also make the bed unstable potentially resulting in injury to the user and patient.

## 6.3 FOLD DOWN SIDE RAILS (OPTIONAL)

Bed rails (optional) are easily fitted by clamping to the bed frame. When fitted correctly, AC4 bed rails meet the requirements of Medical equipment international standard IEC 60601-2-52, which includes specific rail requirements.

### 6.3.1 FITTING THE SIDE RAILS

NOTE: The rails are a left and right hand set and when fitted correctly fold down towards the foot end of the bed.

1. Before fitting the rails unscrew the hand wheel to bring the locking clamp into an open neutral position which will allow the rail to fit over the side rail of the bed.

Closed Position



Open Position



2. With both hand wheels now unscrewed, lift the rail into position. It is important that the rails are fitted to a specific position on the bed, which is indicated by the bolt on top of the bed frame. Position the rail over the bolt and tighten the hand wheel to the closed position.



### 6.3.2 SIDE RAIL OPERATION

The rails are either fully raised, or folded down closed. When the side rail is raised to its full height position it will automatically lock in the positioning hole on the rail clamp plate. The release/lock knob feature can be both visually and audibly confirmed.


**WARNING:** Do not operate side rails when there is potential of entrapment or crushing to the patient, caregivers, visitors and/or other equipment.

#### Side Rail Operation


To raise the side: Lift the top rail until the locating pin/knob automatically clicks into its locked position at maximum height.

**To lower the side:** Hold the top rail near the foot end of the bed (1) and gently pull out the red operating knob (2) lowering the side towards the foot of the bed whilst controlling the descent speed. As the side begins to be lowered the operating knob can be released.

**WARNING:** Care should be taken whilst operating the side rails to ensure hands and fingers are kept away from potential squeeze hazards. Side rails should only be used after clinical assessment and within organisational policy and guidelines.

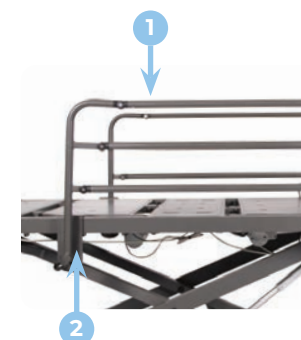

Unlocked

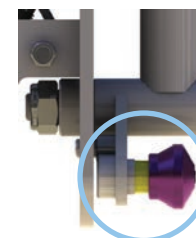

Locked



### 6.3.3 SIDE RAIL SAFETY LOCK INDICATOR

The side rail locating pin has a yellow coating that allows for visual inspection of the locking device. When the side rail is in the raised position, viewing the release knob from the side there should be no yellow shaft visible. If there is yellow shaft visible the safety side is NOT securely locked in the upright position and needs to be adjusted.

The yellow coating should NOT be visible from above when the safety side is correctly locked in position.

*Side rail release and lock knob shown in unlocked and locked positions*

**WARNING:** After adjustment of the side rail visually check the locking pin location. There should be no yellow shaft on the locking pin visible if locked correctly. DO NOT use other manufacturer's side rails on AC4 beds.




## 6.4 KING SINGLE EXTENDABLE FOOT-END MATTRESS RETAINER

AC4 beds come standard with an extendable foot-end mattress-retainer that is used in conjunction with the bed extension.



To extend the mattress retainer hold vertically and insert ends into 1 of the 3 receiving keyhole mounts which are 100mm spaced at the foot end of the mattress platform to match the desired bed extension (nil, 100mm, 200mm). Lower the retainer to the mattress platform. It will now be locked in place at the mount keyholes.

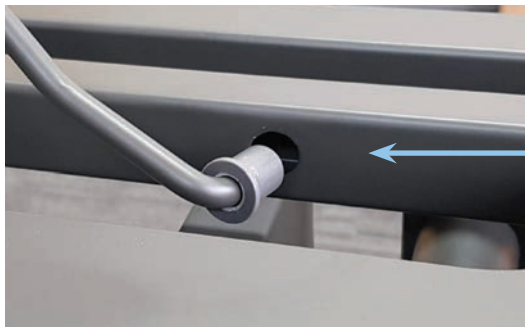

Retaining Pin



## 6.5 BATTERY BACKUP INSTALLATION INSTRUCTIONS

### You will need:

- Phillip head screwdriver
- Flat object such as flat head screw driver or large coin

### Instructions:

1. Battery backup is installed on plate next to control box. Secure the battery by screwing the four screws through the battery into the pre-drilled holes on the plate.

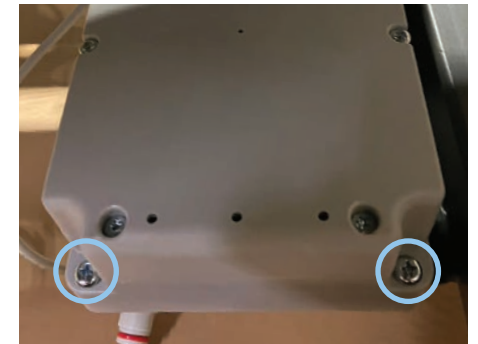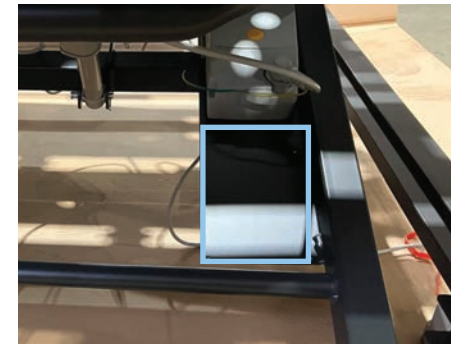

2. Use your flat object to unlock the control box case as shown.




3. Locate the correct port for the battery backup and remove the rubber plug.



4. Plug the battery backup into the port and replace the cover.



## 7. SPECIFICATIONS

SIZE	SINGLE	KING SINGLE
BASIC DIMENSIONS		
Overall Length	2100–2300mm	
Overall Width	920mm	1070mm
Overall Width (rails fitted)	1030mm	1200mm
Mattress Surface Length (between panels)	2000–2200mm	
Bed Platform Height Range	300–750mm	
Back Rest Angle Adjustment	0° to 70°	
Knee Bend Angle Adjustment	0° to 40°	
MATTRESS SIZE		
Recommended Mattress Size Maximum Mattress Depth	Width = 900mm	Width = 1050mm
	Length = 2000mm	
	Depth = 150mm	
BED WEIGHT & SAFE WORKING LOAD		
Without Side Rails	112 kg	116 kg
SWL	220 kg (36.4st)	
ELECTRIC PARAMETER		
Power Input	100–240V AC 50Hz, Max current 1A	
Rate Load	Interruption 10%, Max 2 minutes / 18 minutes	
Leakage Protection	Class 1 Type B 	
Waterproof Grade	IP54	

### STANDARD

Aidacare AC4 beds are designed and manufactured in compliance with International Standard IEC 60601-2-52.

## 8. CLEANING

### WARNING

Before cleaning always disconnect the bed from the mains power source.

Note: It is strongly recommended that protective clothing be used when cleaning.

Note: The following procedure also applies to the beds accessories, but does not include mattresses.

- Wipe the beds surfaces with a soft cloth moistened with hot water and mild detergent (or the facilities recommended cleaning solution). It is recommended that the cloth be a white material to prevent bleeding of dye colours to the bed. Take extra care in areas that can harbour dirt or dust.
- Rinse with clean water and a dry clean cloth.
- To remove potentially infectious materials such as body fluids, or when the bed has been used by someone with a known infection. Clean with NaDCC (e.g. Presept, Actichlor) at 10,000 ppm of available chlorine. Rinse and dry.
- Allow the cleaned parts to dry before replacing the mattress.

### CAUTION

Do not use abrasive Hypocarbonate or Phenolic compounds or pads.

## 9. MAINTENANCE

To prolong the working life of the bed and for trouble free operation it is recommended that routine inspections should be performed by a person who is suitably qualified and well acquainted with the design, use and care of AC4 beds.

If inspection reveals that the safety of the bed is jeopardised in any way, it should be IMMEDIATELY WITHDRAWN FROM SERVICE, tagged as such and the owner notified.

The following preventative maintenance procedures should be carried out every 12 months.

Electrical test & tag – where lifters are in use in a commercial environment then the testing regime applied should align with the risk assessment made by the owner for safety and other compliance requirements.

**NOTE:** Other than general maintenance activities all other service or repair work must be performed by qualified and trained persons who are approved by Aidacare.

**WARNING:** The bed must be disconnected from the mains power supply before starting any maintenance activity.

General	
1	Examine all the warning signs on the bed, examine the bed for obvious signs of damage. All aspects of the equipment should operate as intended. Check that all nuts, bolts and other fasteners are tight and are not missing.
2	Examine flexible cables and conduits for cuts, cracks, abrasions or other deterioration.
3	Check that the power supply plug is not damaged. If either the power cable or plug is damaged then both the cable and plug must be replaced as a complete assembly by an approved service agent.
Central Locking System	
1	Apply the brakes and push the bed forwards. If any of the four braking castors rotate the brake is not fully effective.
2	Check the two steering lock castors for correct operation.
3	If either the brakes or the steering locks are not working correctly contact an approved service agent.

Side Rails (if fitted)	
1	<p>The locking pin mechanism and striker plate arrangement should be checked every three months.</p> <p>Inspect the side rails for general wear and tear on fittings and pivot joints every six months.</p> <p>Check for rust, structural integrity of welds and paint finish every twelve months.</p>
2	<p>Raise the side rail to the upright position until the locking mechanism clicks into place. Refer to the operating procedures starting on page 10 of this manual for further instruction.</p>
3	<p>Check that there is no yellow shaft visible. In the locked position the red on the locking pin shaft should not be visible from above</p>
4	<p>If the red shaft is visible when the side rail is in the raised and locked position, contact Aidacare for service advice.</p>
5	<p>Check that the rails pivot axles are securely fitted to the bed. The pivot axles are locked into place by the manufacturer by a 10mm lock nut.</p>

## 10. TRANSPORT & STORAGE

### Handle With Care.

AC4 beds should be stored in a clean, dry, well ventilated area.

Do not drop and avoid shock or violent impact when transporting.

**The following limits apply during transport and or a storage period of up to 15 weeks duration:**

Ambient temperature	-25°C – 70°C
Relative humidity	10% – 75%
Air pressure	50 kPa – 106 kPa

**The following limits apply to normal operating conditions and or periods of storage longer than 15 weeks:**

Ambient temperature	10°C – 40°C
Relative humidity	30% – 75%
Air pressure	70 kPa – 106 kPa

## 11. ACCESSORIES

A range of AC4 accessories are available. In addition to the below listed accessories, Aidacare specializes in other complimentary products such as foam mattresses, pressure care surfaces, over bed tables, side lockers and adjustable day chairs.

Please contact Customer Service on **1300 133 120** for further information.

Product code	Description
BEA829100	3/4 Folding Side Rails
BEA829105	Clamp On Bed Rail
BEA829140	Self Help Pole
BEA029710	IV Pole
BEA829150	Padded Surround (Single)
BEA829160	Padded Surround (King Single)
GSP405130	Battery Backup



Do not use accessories that have not been designed for use with the AC4 beds.

## 12. WARRANTY & SERVICE

Aidacare AC4 Beds carry a comprehensive warranty as detailed below. Please keep proof of purchase for proof of warranty commencement.

**10** YEARS WARRANTY  
Frame

**3** YEARS WARRANTY  
Electronics

**1** YEARS WARRANTY  
Handsets & Cable

### Consumer Warranty Against Defects:

Aidacare is bound by Australian Consumer Law in respect to Warranty Against Defects. A copy of the Aidacare Warranty Against Defects statement is available on the **Aidacare homepage** at **aidacare.com.au**. All warranty matters should be directed to **Aidacare Customer Service** on **1300 133 120**.

Aidacare does not warrant against excessive or incorrect use, modification or any situation that could not be deemed as fair wear and tear. This is a back to base warranty and does not cover freight costs pertaining to the return of any items under warranty. Aidacare acknowledges and conforms to all ACCC guidelines around consumer rights. For further information contact **Aidacare** on **1300 133 120** or email **product@aidacare.com.au**.

Aidacare will not warrant the safety and or correct functioning of products where any original components have been changed or modified by non-Aidacare approved and trained service & maintenance staff or external providers. Furthermore, safety is not guaranteed where components have been replaced with non-original Aidacare approved parts. If any faults are detected upon receipt of this product please phone Aidacare. Any faults that are detected during normal use should be reported to Aidacare immediately to determine if warranty conditions apply and if so, the necessary repair or replacement work to be completed.

Spare parts list is available upon request.

For servicing, preventative maintenance and any other questions regarding this or any other Aidacare product please contact:

### AIDACARE PTY LTD

Building 3A, 1 Moorebank Avenue,  
Moorebank NSW 2170  
Australia

**1300 133 120** | [aidacare.com.au](http://aidacare.com.au)

### AIDACARE LIMITED

Arcadia House, Maritime Walk,  
Ocean Village, Southampton SO14 3TL  
United Kingdom

**01622 541235** | [aidacare.co.uk](http://aidacare.co.uk)

## This image shows a single sheet of white paper with horizontal blue ruling lines. The lines are evenly spaced and run across the width of the page. There is no text or other markings on the paper.

## **AIDACARE PTY LTD**

Building 3A, 1 Moorebank Avenue,  
Moorebank NSW 2170  
Australia

**1300 133 120** | [aidacare.com.au](https://aidacare.com.au)

## **AIDACARE LIMITED**

Arcadia House, Maritime Walk,  
Ocean Village, Southampton SO14 3TL  
United Kingdom

**01622 541235** | [aidacare.co.uk](https://aidacare.co.uk)

MADE IN PRC