

MATTRESS RANGE

Benefits

Australian made covers using the highest quality European two way stretch fabric.

- Australian made foam
- Fire retardant

Removable Cover:

1. Easy inspection
2. Facilitates better hygiene
3. Fully washable


Mattresses do not require rotation, minimising OH&S handling risks.



LIFECOMFORT

Dual layer support surface incorporating a full size memory foam overlay and durable high density foam base. European fabric cover has fully welded seams and water resistant zippers with reinforced waterfall flaps for maximum infection control.


Code	Description	Length	Width	Height	Warranty	S.W.L
BEM046410	Lifecomfort - AC Series Single	2000mm	900mm	130mm	5 Years Foam 2 Years Cover	120kg
BEM046425	Lifecomfort - Small	1920mm	820mm	130mm	5 Years Foam 2 Years Cover	120kg
BEM046460	Lifecomfort - King Single	2030mm	1050mm	130mm	5 Years Foam 2 Years Cover	120kg

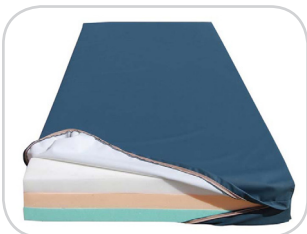

LIFECOMFORT DELUXE

Three layer PUPPS compliant support surface incorporating a full size memory foam overlay and two layers of high quality foam for maximum comfort and support. European fabric cover has fully welded seams and water resistant zippers with reinforced waterfall flaps for optimum infection control.

Code	Description	Length	Width	Height	Warranty	S.W.L
BEM046510	Lifecomfort Deluxe - AC Series Single	2000mm	900mm	150mm	5 Years Foam 2 Years Cover	150kg
BEM046525	Lifecomfort Deluxe - Small	1920mm	820mm	150mm	5 Years Foam 2 Years Cover	150kg
BEM046455	Lifecomfort Deluxe - King Single	2030mm	1050mm	150mm	5 Years Foam 2 Years Cover	150kg

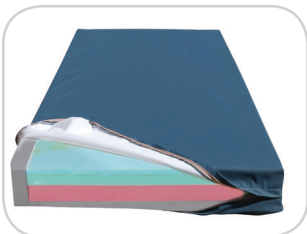

LIFECOMFORT BARIATRIC

Four layer foam, designed for larger users with a maximum user weight of 300kg. Full size memory foam overlay and side walls for added stability and ease of transfers. European fabric cover with fully welded seams, water resistant zippers and reinforced waterfall flaps for optimum infection control.

Code	Description	Length	Width	Height	Warranty	S.W.L
BEM046700	Lifecomfort Bariatric - King Single	2030mm	1050mm	200mm	5 Years Foam 2 Years Cover	300kg



LIFECOMFORT PRESSURE REDUCING

Top of the range PUPPS2 compliant mattress featuring a quad layer support surface with additional side wall supports. European fabric cover has fully welded seams and water resistant zippers with reinforced waterfall flaps for optimum infection control.

Code	Description	Length	Width	Height	Warranty	S.W.L
BEM046610	Lifecomfort Press Red - AC Series Single	2000mm	900mm	150mm	5 Years Foam 2 Years Cover	150kg
BEM046625	Lifecomfort Press Red - Small	1920mm	820mm	150mm	5 Years Foam 2 Years Cover	150kg
BEM046660	Lifecomfort Press Red - King Single	2030mm	1050mm	150mm	5 Years Foam 2 Years Cover	150kg
BEM046675	Lifecomfort Press Red - Double	2030mm	1370mm	150mm	5 Years Foam 2 Years Cover	150kg
BEM046680	Lifecomfort Press Red - Queen	2030mm	1520mm	150mm	5 Years Foam 2 Years Cover	150kg

MATTRESS RANGE



ASSET MANAGEMENT

Serial numbers, heat transferred onto the cover's base material, allow tracking on fixed asset registers.



INFECTION CONTROL


1. Fully welded seams
2. Reinforced waterfall flaps
3. Water resistant zippers



LIFECOMFORT WITH SAFETY SIDES

This mattress features unique side walls that extend up above the height of the mattress providing security for the patient. European fabric cover has fully welded seams and water resistant zippers with reinforced waterfall flaps for optimum infection control.

Code	Description	Length	Width	Height	Safety Sides Height	Warranty	S.W.L
BEM046475	Lifecomfort with Safety Sides	1940mm	900mm	130mm	220mm	5 Years Foam 2 Years Cover	120kg

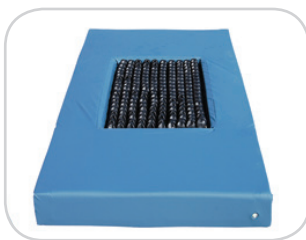

LIFECOMFORT PRESSURE REDUCING WITH ROHO CUT OUT

Designed to accept the standard size ROHO mattress section (BEM636320). This allows for additional pressure relief in the pelvic section of the mattress providing the ultimate in support, comfort and pressure reduction all in one. The foam mattress also features a 4 layer foam design, side walls and two way fully welded waterproof stretch cover.

Code	Description	Length	Width	Height	Warranty	S.W.L
BEM046720	Lifecomfort Pressure Reducing with ROHO Insert Cut Out	1940mm	900mm	150mm	5 Years Foam 2 Years Cover	150kg

 Lifecomfort Mattress Range is registered in the ARTG. AUST R ARTG 288486


HIGH QUALITY MATERIALS

The Lifecomfort Mattress range is manufactured in Australia using the best quality materials available.



FOAM:

The foams used in our mattresses have passed the most stringent industry testing including Pressure Mapping, using the X-sensor device. All of our mattresses are treated with Purify™, a performance anti-microbial which is incorporated into foam during manufacture to produce a long-lasting, cleaner and fresher product. Purify™ helps prevent and control odour-causing bacteria and degradation that can be caused by mould and fungi.

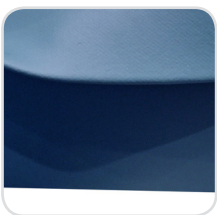

COVERS:

"CARRFLEX Bi-elastic" is a PU coated fabric having a fire retardant and antibacterial polyurethane coating on a 2-way stretch polyamide base fabric. It has been specifically designed for use in the medical mattress market.



Australian Made