





# REHAB RX JNR TILT-IN-SPACE WHEELCHAIR

Designed to provide advanced postural support and pressure care for children. It offers various adjustable features for improved positioning and comfort, making it suitable for children with specific care needs.

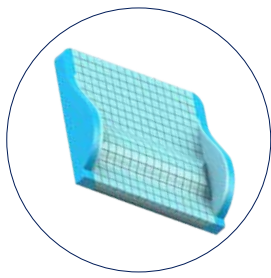

**CONTOURED  
BACKREST**

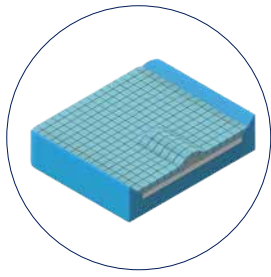

**CONTOURED  
SEAT**

CRASH TEST  
**APPROVED**  
AS/NZS 3696.19

**1** YEAR WARRANTY  
Components

**2** YEARS WARRANTY  
Frame

SWL  
**180kg**

SWL  
**23.8st**



### Sturdy Frame

Rigid steel construction for maximum durability

### Quick-Release Wheels

Solid 20" tires with easy height and depth adjustment for versatile terrain handling

### Height-Adjustable Castors

6" wheels for smooth, adaptable maneuverability

### Dual Braking System

Attendant-operated drum brakes and push-to-lock handbrakes for enhanced safety

### Ergonomic Headrest

Height, depth, angle, and width adjustments; fully removable for user flexibility

### Breathable Upholstery

PU multi-stretch fabric covers for comfort and airflow

### Premium Cushioning

Multilayer castellated back and seat cushions in breathable PU for pressure relief

### Lightweight Design

Starts at 39.5 kg for easier transport and handling




## ADJUSTMENTS




Backrest Height



Seat & Backrest Depth

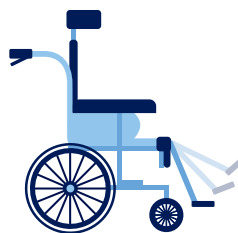

Elevating Legrest

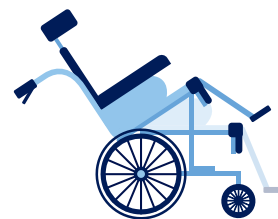

Tilt-In-Space  
26° Tilt | 30° Recline

## Discover Our Exclusive Range of Accessories

IV Pole  
WCA848280

Safety Lap Belt  
WCA848300

Oxygen Bottle Holder  
WCA848290

Tray Table  
WCA848570

Stump Support  
Various sizes

REHAB RX JNR WHEELCHAIR	Seat Width	Seat Depth	Front Seat Height	Rear Seat Height	Overall Width	Backrest Height	Armrest Height	Min. Chair Weight	SWL
	mm / in	mm / in	mm / in	mm / in	mm / in	mm / in	mm / in	kg / lbs	kg / st
MWS449630	320 - 370/ 13 - 15	440 - 490/ 17.3 - 19.29	410 - 440/ 16.1 - 17.3	400 - 475/ 15.8 - 18.7	600 / 23.6	460 - 710 / 18 - 28	310 - 460/ 12.2 - 18.1	39.5/ 87	180/ 28.3

\*Range of motion dimensions are accurate to a +/-10mm tolerance

ASPY142.1 241024