



# ULTIMATE 8 MATTRESS REPLACEMENT

## SUITABLE FOR:

- Aged Care
- Home Care
- Hospital




The Ultimate 8 mattress caters to users considered to be at high risk of pressure injury. Unique pump technology reduces noise, promoting improved sleep quality. The mattress is easy to clean and designed for durability.



Infection Control



Skin Microclimate

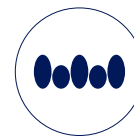

Alternating Pressure Relief



Easy Maintenance & Cleaning



Pressure Redistribution



Adjustable User Weight Settings

### HIGH LEVEL PRESSURE THERAPY:


- Promotes circulation and provides pressure redistribution to treat and prevent pressure injuries
- Rapid release CPR allows quick deflation of the mattress, including head cells, in emergencies
- Prevents bottoming out by increasing sacral zone air pressure when the mattress backrest is raised
- 24hr transport mode and static support layer create a security buffer allowing patients to remain on an inflated surface in the event of power failure
- Easy-to-read LCD screen provides clear and precise information on cycle time, operating mode, and fault codes to assist with troubleshooting
- Improves workplace safety via integrated power cord management in the mattress cover
- Designed for durability, the cell guard system improves mattress strength in key entry/exit areas

### IMPROVES SLEEP QUALITY:

- Provides a whisper-quiet sleep environment using innovative pump to reduce noise and vibration
- Improves sleep comfort by passing room temperature air beneath the sleep surface, helping to regulate body temperature

### EASY INFECTION CONTROL:

- Prevents fluid and bacterial ingress with welded seams and zip protection
- Multi-stretch cover reduces shear and offers a high moisture vapour transfer rate (MVTR) to control skin microclimate
- Easy to clean mattress cover can be laundered at up to 85°

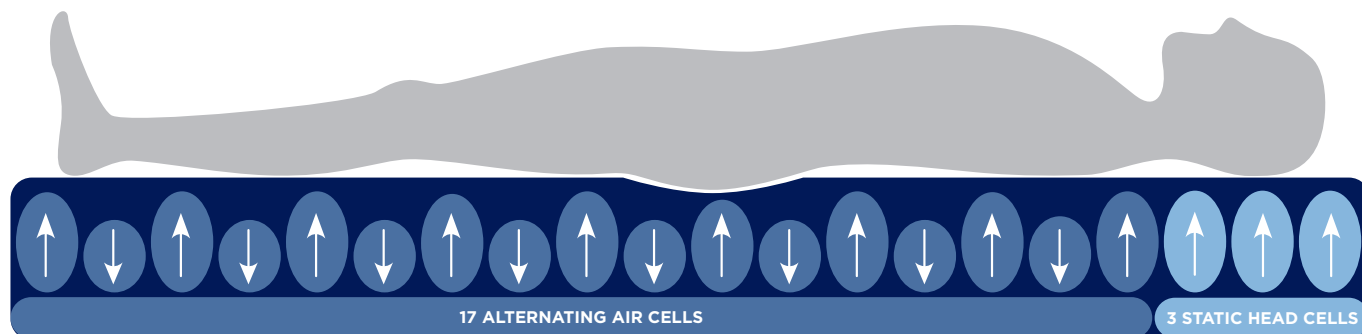

Pressure redistribution during alternation cycle



ULTIMATE 8 MATTRESS	SIZE (mm)	OVERALL LENGTH (mm)	OVERALL WIDTH (mm)	OVERALL HEIGHT (mm)	NUMBER OF CELLS	STATIC HEAD CELLS	CELL CYCLE	CYCLE TIME (mins)
PTM639320	Single	2000	900	200	20	3	1 in 2	11
PTM639330	King Single	2000	1050	200	20	3	1 in 2	11

ASFY85.11\_080822