

NEW PRODUCT LAUNCH




Whill C2 Custom Backrest Cane Assembly - Now Available!

We have designed a new adjustable backrest cane that allows customers to fit any aftermarket backrest cushions such as Jay, SPEX, or Varilite.

These canes are simple to install and can be width adjusted to accommodate backrests from 16" to 20" in width.

Standard Configuration is 300mm in height and allows the installation of most common backrest heights. If there is a need to install a taller backrest, our design will enable you to simply swap the vertical upright with a 500mm tube to increase the height of the backrest cane assembly.

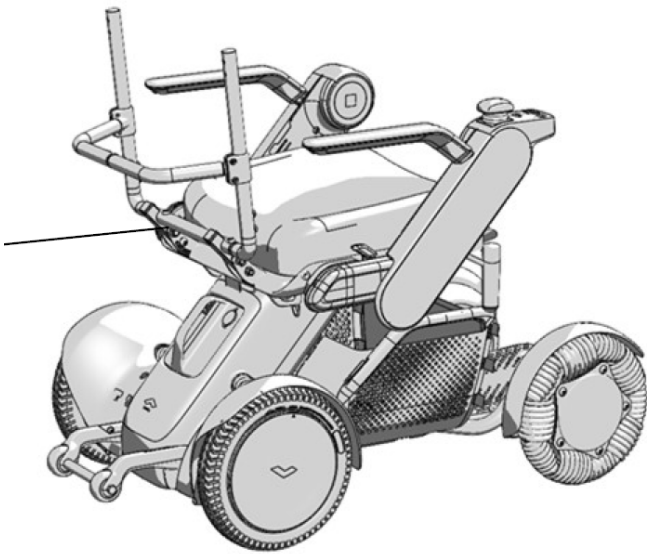

Custom backrest Canes – Ready to install aftermarket Backrest cushions



Assembled with Custom Backrest Cushion

NEW PRODUCT LAUNCH

PRODUCT CODES

PRODUCT FEATURES

- Install any aftermarket backrest cushions - Allows users to use any complex backrest cushions
- Width Adjustable – Install backrest cushions from 16” to 20” width
- Height adjustable – Taller backrest upright allows installation of taller custom backrest cushions

INSTALLATION & SUPPLY

- Please refer to the video below for detailed fitment instructions for the custom backrest cane assembly onto the Whill C2 Wheelchair (written instructions are also supplied with stock)
 - [Model C2 Adjustable Backrest Mounting System \(youtube.com\)](#)
- For mounting the backrest itself, as each brand differs, salespeople are encouraged to work with their local CASE teams
- These cane systems do not include the cost of an aftermarket backrest, this must be quoted and sold separately
- Where users require the taller cane option, both WCA870456 & WCA870457 need to be added to the sales order/quote